

Controlled precision




FLUID CONTROL Arthro

Irrigation and Suction pump for arthroscopy
and full-endoscopic spine surgery

Performance in a new dimension

FLUID CONTROL Arthro is a new pump system for arthroscopic endoscopy and full-endoscopic spine surgery. The particular strength of this system lies in a partial combination of high performance and cost-effective overall concept. The new software and technical enhancements make FLUID CONTROL Arthro ideal for diagnostic and therapeutic interventions.

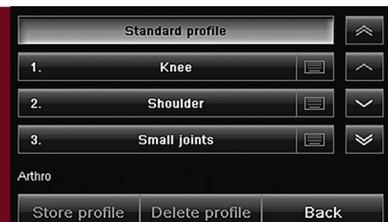

Reserve capacity for all cases

The system is equipped with high pressure and pumping capacities, and this allows it to offer adequate capacity reserves for all eventualities. This enables optimum visibility to be quickly reinstated specifically in wash mode.

- Pressure: up to 200 mmHg
- Pumping capacity: up to 2 l/min

Intelligent fluid management

The particularly refined control quickly equalizes the pressure but also recognizes abrupt interruptions in the irrigation flow, for example by closing the stopcock on the endoscope. Optical and acoustic alerts are given when critical overpressures occur.

Particularly user friendly

It is not necessary for application parameters to be continually re-entered as a result of the pre-defined and freely programmable indication profiles.

Spine Mode

The spectacularly safe flow control was specially developed for use in full-endoscopic spine surgery:


- Optimum endoscopic visualization of the spinal canal and disk space
- Improved safety through continuous pressure monitoring and over-pressure control
- Standardized application profiles for VERTEBRIS instrument systems
- Low consumption of irrigation fluid through optimally regulated flow performance

Innovative tubing concept

Tubing systems from Richard Wolf give you an optimum combination of user comfort, safety and cost effectiveness.

Complete tube set

Once a day or less frequent use: we recommend a **complete tube set**. Make your choice here between the re-usable and single-use version.

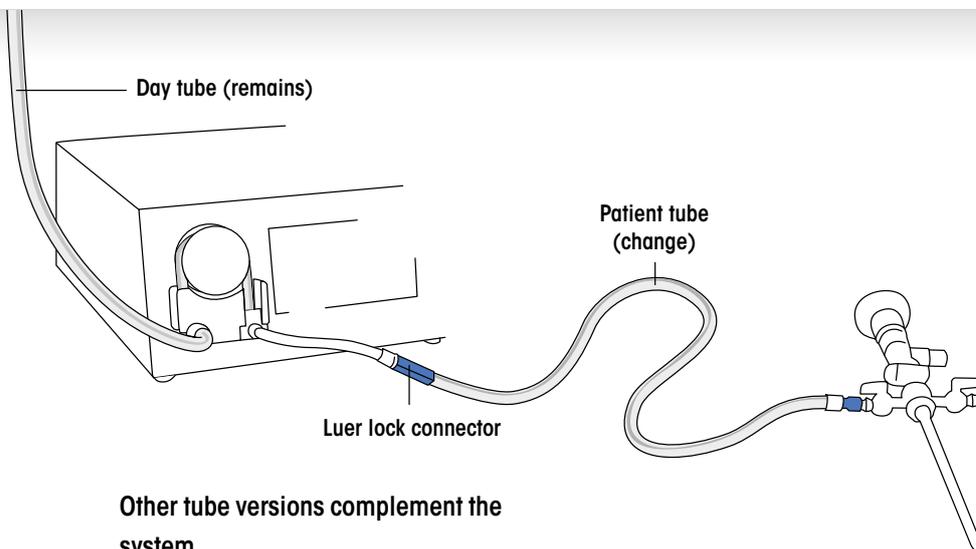

Other tube versions complement the system.

Day-patient tube set

Already from two interventions/patients each day: make use of the **day-patient tube set**. This set is made up of a day tube which remains in the pump throughout the day and a patient tube which is changed for each intervention. The irrigation fluid bad can remain on the pump until it has been completely used up. **This significantly reduces the amount of waste and cuts costs.**

Note: The use of the day-patient tube set has not been approved for the application spine mode.

Safety is our commitment

FLUID CONTROL Arthro keeps track of components. The system is fitted with intelligent **transponder technology**. This automatically recognizes whether the tube is a system-compatible Richard Wolf tube, the tube type and the number of applications still permitted.

- Reduces the risks of infection
- Optimizes the comfort factor for personnel

Single genius – Ingeniously doubled

FLUID CONTROL Arthro **combines the irrigation and suction pump** in a single unit. A vacuum pump is integrated alongside the roller pump for irrigation. This makes significant cost savings.

- Two pumps in one unit
- A single tube is adequate for evacuation. A second roller-pump tube set is not required.

Comprehensive user comfort


The new equipment design with sloping front panel makes the large touchscreen conveniently accessible.


- Intuitive operation
- Autoclavable remote control
- Foot switch for activation of wash mode

FLUID CONTROL Arthro has already been prepared for **integration in core**, the digital Operating Room from Richard Wolf.

Ordering data

Pump and accessories	
	Set Arthroscopy: Arthro/Spine suc./Irrig. Pump 200 mmHg (2204001), irrigation tube set spike, L 3 m (4171223), power cable (2440.03)
	Set Spine Surgery/Arthroscopy: Arthro/Spine suc./Irrig. Pump 200 mmHg (2204001), spine mode software module (2204101), irrigation tube set spike, L 3 m (4171223), power cable (2440.03)
	Footswitch2204901*

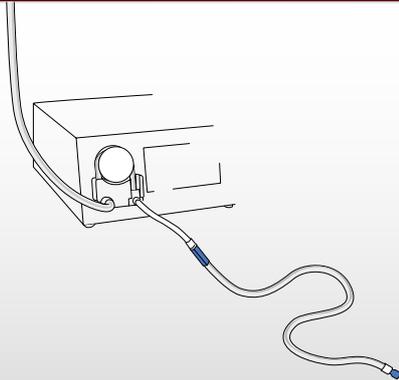Irrigation	
	Accessories - reusable Irrigation tube set spike , incl. 10 replacement sealing membranes, silicone, PACK=1 PC, with 2 piercing spikes, 20x reusable, for Fluid Control Arthro/Spine 2204 and Fluid Control Lap 2216, reusable 8171223
	Accessories - for single use Irrigation tube set spike L 3 m , with Y-piece, PVC, PACK=10 PCS, for Fluid Control Arthro/Spine 2204, Fluid Control Lap 2216 and Fluid Control 2225, sterile, for single use 4171223
	Irrigation tube set care lock L 3 m , with Y-piece, PVC, PACK=10 PCS, for Fluid Control Arthro/Spine 2204, Fluid Control Lap 2216 and Fluid Control 2225, sterile, for single use 4171224
	Day irrigation tube set , with Y-piece, PVC, PACK=10 PCS, for Fluid Control Arthro/Spine 2204, sterile 4171226*
	Patient irrigation tube set L 3 m , PVC, PACK 10x patient irrigation tubing set, with Luer Lock, for Fluid Control Arthro/Spine 2204, sterile 4171227*

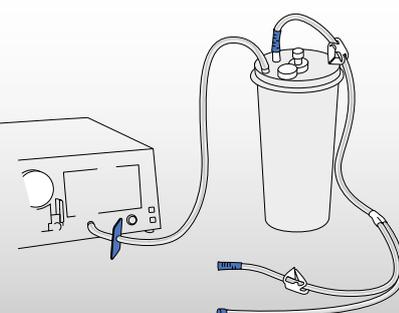Suction / Evacuation	
	Accessories - reusable Tube set for vacuum, incl. filter , useful life of 30 days, unsterile 2206207
	Secretion container 2 L , with cover, reusable 8170655
	Accessories - for single use Tube set for suction (with 2 ports) , single-use, sterile, (PACK=10 PC.), (with Y-piece, for connecting the suction container with shaver handle/trocar sleeve) T0503-01

* Not approved for application spine mode